


2.

Upgrade existing **mobile** X-ray systems

Mode of operation using Toughpad/Tablet-PCs

Wireless DR flat panel

Automatic synchronisation of detector and generator by means of **AED***
No need to modify the X-ray system and to adjust the system or the cable connections

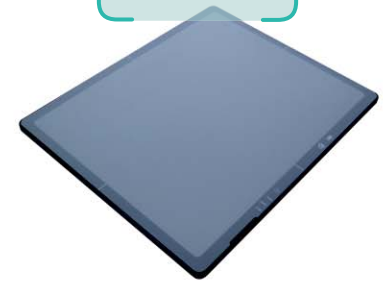

Mobile X-ray unit with **operating console**




Transfer of raw images to Toughpad or Tablet-PC – no additional access point required

Operating console Toughpad or Tablet-PC with and without keyboard with worklist (DICOM worklist)



A variety of mounting systems are available for the operating console

Transfer of images to the practice network via wireless LAN

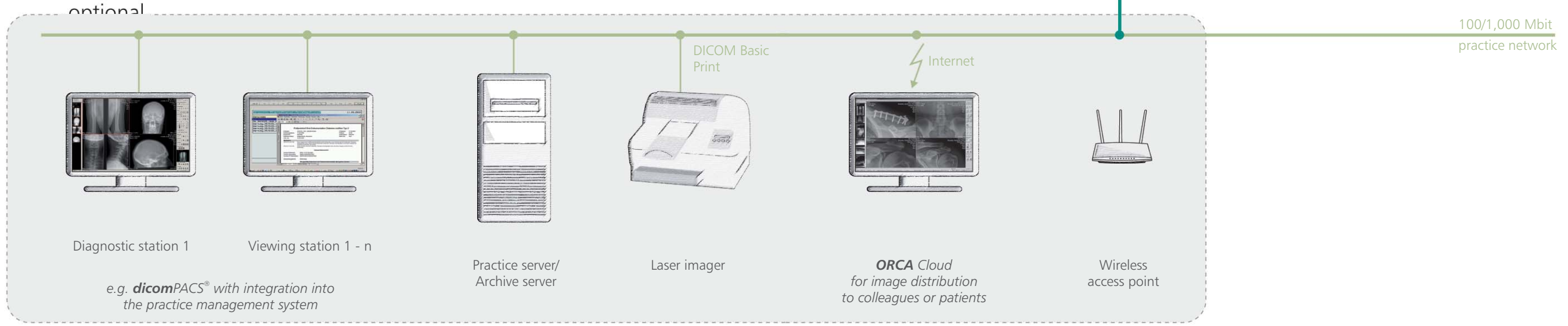

* Auto Exposure Detection